

1941400000000/01.22/v002
© 2022 – UVEX ARBEITSSCHUTZ GMBH

UVEX ARBEITSSCHUTZ GMBH
Wuerzburger Str. 181-189
90766 Fuerth
GERMANY
T: +49 911 9736-0
F: +49 911 9736-1760
E: safety@uvex.de
I: uvex-safety.com

a Marking on frame

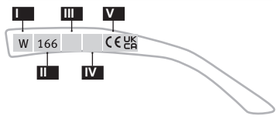b Marking on lens

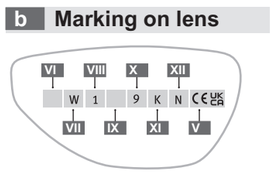uvex

uvex safety eyewear

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UKCA
Approved Body Number 2571
TÜV Rheinland UK Ltd
Friars Gate (Third Floor), 1011 Stratford Road, Shirley, Solihull, B90 4BN
UNITED KINGDOM

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EN Instruction leaflet for uvex safety eyewear according to EN166:2001, EN167:2001, EN168:2001, EN169:2002, EN170:2002, EN171:2002, EN172:2001, GS-ET29:2019, ANSI/ISEA Z87.1-2020 (USA). Certain models comply with CSA Z94.3. Please ask a uvex specialist which products are included.

This safety eyewear satisfies the requirements of PPE Regulation (EU) 2016/425 and Directive 2001/95/EC and the Personal Protective Equipment Regulations (Regulation (EU) 2016/425/ as brought into UK law and amended). The Declaration of Conformity can be found at:
EU: uvex-safety.com/ce **UK**: www.uvex-safety.co.uk/en/media/uk-declarations-of-conformity
The model name or number is indicated on the label or on the safety eyewear itself. Transport and store in the original packaging only, in a dark, dry location where the product will not be subjected to mechanical loads, at a temperature not exceeding 40°C and at <70% air humidity. Date of manufacture (MM/YYYY): on the safety eyewear and packaging. Shelf life if stored correctly: six years from the date of manufacture. Recommended service life: max. two years from first use, depending on the intensity of use. Before each use, check the safety eyewear for damage and ensure that it fits properly. The safety eyewear only provides protection for the area that it covers. If necessary, you may need to wear sealed safety goggles or a face shield. Do not use but replace instead eye and face protection that have been subjected to impact, is scratched, worn or damaged. Clean the eyewear as regularly as possible using uvex cleaning fluid or by rinsing it under running water (do not clean while it is dry); clean immediately if it comes into contact with chemicals. When using commercially available disinfectants, only spray or wipe down the eyewear. Use only original uvex spare parts and accessories. Persons with sensitive skin may suffer an allergic reaction upon contact with certain materials. If the safety eyewear is worn over prescription glasses, mechanical force can be exerted on the glasses and the wearer, potentially causing damage or a hazard. Mechanical protection at extreme temperatures (−5°C/+55°C) is indicated with the additional T certification mark. If this marker is not present, the safety eyewear is only suitable for protecting the wearer against high-speed particles at room temperature. If the lenses and the frame do not bear the same code (S, F, B, A), the glasses should only be used in environments corresponding to the lower category. May catch fire in case of direct contact with open flames or hot surfaces. To protect against optical radiation, filters that are appropriate for the application must be used. These are labelled with a prefix followed by the protection level (VI).

- EN169:2002 (welding protection filter without prefix – level of protection only)
- EN170:2002 (UV protection filter, prefix 2 or 3)
- EN171:2002 (IR protection filter, prefix 4)
- EN172:2001 (sun protection filter, prefix 5 or 6)

Recommendations on the use of filter categories can be found in the relevant standards or the uvex catalogue. Some applications require protection against mechanical, chemical, thermal or electrical hazards. So that the correct eyewear is selected, a risk assessment must always be carried out by the health and safety manager. Do not use filters with a transmittance of <75% at dusk or during the night. Only filters with signal light detection (additional code C or in accordance with EN172:2001) are suitable for use in traffic. This eyewear is not suitable for use as laser protection.

a	Marking on frame
I	Manufacturer identification
II	Number of the EN standard
III	Field(s) of use
IV	Mechanical strength
V	Certification mark

IV	IX	Mechanical strength
none		without mechanical strength
S		increased strength
F		low energy impact (45 m/s)
B		medium energy impact (120 m/s)
A		high energy impact (190 m/s)
T		extreme temperatures (−5 °C/+55 °C)

b	Marking on lens
VI	Level of protection
VII	Manufacturer identification
VIII	Optical class
IX	Mechanical strength
X	Non adherence of molten metal and resistance to penetration of hot solids
XI	Resistance to surface damage by fine particles
XII	Resistance to fogging
V	Certification mark

III	Field(s) of use (where applicable)
none	General use, non-specific mechanical risks. Risks due to UV, IR and/or visible radiation
3	Liquids (droplets and splashes)
4	Coarse dust particles: Grain size >5 µm
5	Gas, vapour, mist, smoke and fine dust particles: Grain size <5 µm
8	Arc flash due to short circuit in electrical equipment
9	Splashes of molten metal and penetration of hot solids

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